

NEVI™

Sickle Scaler Collection



Instrumentation and
Procedural Products



Innovative Designs
For Effective Deposit
Removal

 **HuFriedyGroup**
The Best In Practice

Nevi™ Sickle Scalers

Better Than Ever: EverEdge™ 2.0 Technology

The Nevi™ Series features 5 scalers capable of scaling light to heavy calculus with minimal tissue trauma and superior adaption to allow precise access in difficult areas.

The Nevi™ Instruments were specially designed with input from clinicians—dental hygienists like you who know: the details of oral anatomy, how the right bend or curve of the working end will allow you to successfully scale particular tooth surfaces, and that using the best instruments is what their patients deserve.

Each Nevi™ Scaler was created from a combination of customer feedback and our tradition of quality, innovation and skill. The unique angles and distinct characteristics of each Nevi™ Instrument allow for access and reach into difficult areas, while our **EverEdge™ 2.0 Technology** keeps the instruments sharper longer—adding up to a more efficient, comfortable scaling experience.

Stays
Sharper
Longer

Resin Eight™
Handle

Harmony™
Ergonomic
Handle

Popular #9
Handle

Anterior Sickle Scalers

NEVI™ 1



NEVI™ 1/H5




Two distinct working ends provide solutions for anterior scaling needs

- **Sickle End:** A moderate length, thin scaler with two cutting edges and a blade that tapers to a point. Ideal for removing anterior interproximal calculus and stain
- **Disc End:** A moderate sized, sharp rounded disc. Excellent for removal of anterior lingual stain

The Nevi 1™ Scaler is also available with a longer sickle working end (Nevi1/H5).



Mesial Left Central Incisor

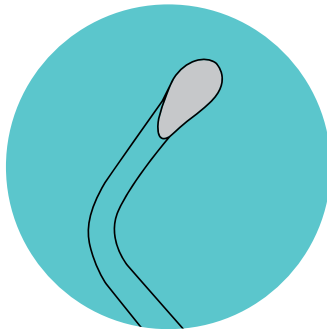

Nevi™ 1 Disc End



Lingual Left Lateral Incisor

Posterior Sickie Scalers

NEVI™ 2



- Perfect for **light to moderate interproximal** calculus removal
- The thinnest dual tapered blades in the Nevi™ Instrument collection remove calculus with minimal tissue trauma
- Deep contra-angles are paired with extra blade length for easy access under contacts



Mesial Right First Molar



Distal Left Second Molar

NEVI™ 3



- Used for light to moderate calculus removal in **tight, narrow posterior interproximal areas**
- Ultra slender, dual cutting blades for effective placement
- Contra-angle bends in shank offer superior reach and adaptation to posterior interproximal areas



Mesial Left First Molar



Distal Left Second Premolar

Posterior Sickle Scalers

NEVI™ 4



- Ideal for scaling **heavy posterior interproximal** deposits
- Increased shank size and rigidity generates effective removal of tenacious, heavy and burnished deposits
- Working ends with dual tapered cutting edges allow for precise interproximal access
- Largest blade in the Nevi™ Instrument collection



Mesial Right First Premolar

NEVI™ 4 Lite




- 10% thinner than the standard Nevi™ 4 Blade for enhanced access to tight interproximal spaces



Distal Right Second Molar

Nevi™ Sickle Scalers Featuring EverEdge™ 2.0 Technology

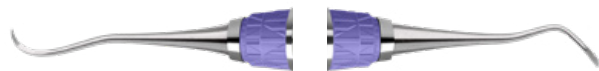

Nevi™ 1 Anterior Sickle Scaler

Harmony **SCNEVI1XE2**

#9 **SCNEVI19E2**

#C8, Lavender **SCNEVI1C8E2**




Nevi™ 1/H5 Anterior Sickle/Discoïd Scaler

Harmony **SNEV1/H5XE2**

#9 **SNEV1/H59E2**

#C8, Lavender **SNEV1/H5C8E2**

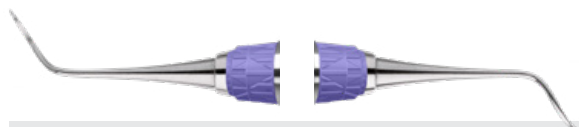


Nevi™ 2 Posterior Sickle Scaler

Harmony **SCNEVI2XE2**

#9 **SCNEVI29E2**

#C8, Red **SCNEVI2C8E2**

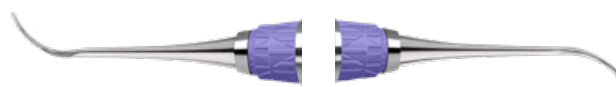

Nevi™ 3 Posterior Sickle Scaler

Harmony **SCNEVI3XE2**

#9 **SCNEVI39E2**

#C8, Red **SCNEVI3C8E2**




Nevi™ 4 Posterior Sickle Scaler

Harmony **SCNEVI4XE2**

#9 **SCNEVI49E2**

#C8, Red **SCNEVI4C8E2**

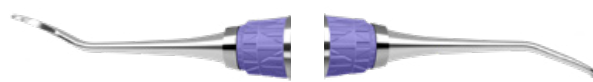

Nevi™ 4 Lite Posterior Sickle Scaler

Harmony **SCNEVI4LXE2**

#9 **SCNEVI4L9E2**

#C8, Red **SCNEVI4LC8E2**

Several Nevi™ Scaler designs are also available with Immunity Steel™ Handles. Please see the full HuFriedyGroup Product Catalog for more information.

To learn more about our Best in Practice scalers and curettes, visit
[HuFriedyGroup.com/Scalers](https://www.hufriedygroup.com/Scalers)

Hu-Friedy Mfg. Co., LLC, 1666 E. Touhy Ave., Des Plaines, IL 60018 | [HuFriedyGroup.com](https://www.hufriedygroup.com)

All company and product names are trademarks of Hu-Friedy Mfg. Co., LLC its affiliates or related companies, unless otherwise noted. Marks not registered in all jurisdictions.

©2024 Hu-Friedy Mfg. Co., LLC. All rights reserved. HF-270/2403

 **HuFriedyGroup**
The Best In Practice